



## Samples preparation guidelines for 10X Genomics 3'GE scRNA-seq v4

*This protocol for the classical 3'GE scRNA-seq serves as our main 10XG scRNA-seq guidelines. Additional guidelines are available for specific applications (VDJ, multiome, multiplexing...), check on our website.*

**CAUTION:** this protocol incorporates changes specific to the [new v4 chemistry](#). In particular:

- Max number of cells that can be targeted is now 20k instead of 10k. Doublets remain low (see table).
- Recommended cell concentration updated from 700-1'200 to 1'000-1'500 cells/ul.
- Increased sensitivity (number of genes detected).
- V4 is compatible with totalseq and "on-chip" multiplexing but not anymore with cellplex multiplexing.

### Disclaimer

*10XG scRNA-seq works in general well but is still technologically at the edge. There is indeed a non-null percentage of wells/runs that fail for purely technical reasons (clogs, instrument issues). Of course, we don't charge for these unless we spot a critical error from the user side (e.g. failure to follow cell prep guidelines/buffers), or if users gave us their green light to carry on despite poor cells quality. Importantly, if a putative failure represents a major issue for your project (e.g. months-long mice generation), tell us in advance and we'll discuss the risks and means to minimize them. As a general rule, if you have extra cells, and if the experiment is critical, keep the cells on ice until you know the instrument run went well, freeze them as a backup with standard techniques, or fix them with DSP/methanol according to 10XG protocol CG000776.*

*For complex/large experiments, we strongly recommend doing a small pilot (1 sample, shallow sequencing).*

## Scheduling the experiment

In general, a 2-3 weeks prior notice is needed for a regular scRNA-seq experiment, but more may be needed if cells isolation must be optimized. In any case, please do not start an experiment that requires samples collection on a specific day or week before having checked/scheduled with us.

Early in the process, send us the submission form to be found on our web site, even if only partially filled in at this stage.

## Interactions with GECF on experiment day

**IN SHORT:** Bring your cells prepared and washed according to these guidelines, on ice, at 1'000-1500 cells/ul, in minimum 50ul, without EDTA, in 1.5-2ml tubes. Bring medium for putative dilutions. Indicate number of targeted cells, and bsl2 level if relevant.

## Targeted cell number

- Indicate the number of cells you want in your data ("targeted" cells number) in the submission form. Recovery rate is uncertain and depends on characteristics of the cells, and on experiment-specific factors, such as too low/too high cell concentration, viability, debris, etc. therefore targeted cells number is only indicative and the number of recovered cells can differ significantly. If you absolutely want a given minimum number of cells, indicate it in the submission form and we will add a safety margin.

- Each additional cell comes with a sequencing cost. Once the cells have been processed, it is not possible to sequence only a fraction of them.
- The ideal number of targeted cells depends on the biological question. For comparing two populations, 2'000 cells may be enough. At the other extreme, for identifying new rare subpopulations of cells (<1%), 15'000 or more cells may be needed. If unsure, we recommend 6'000 cells as a good starting point.
- The doublets rate increases with targeted cells number (see table below), therefore we recommend not targeting more than 12'000 cells unless necessary.

Nb of Recovered Cells	Multiplet Rate (%)
1 000	~0.4%
2 000	~0.8 %
4 000	~1.6%
6 000	~2.4%
8 000	~3.2%
10 000	~4%
12 000	~4.8%
14 000	~5.6%
16 000	~6.4%
18 000	~7.2%
20 000	~8%

## Cells preparation

Cells must be submitted ready for subsequent steps in the procedure, according to the guidelines below and to the latest version of 10XG Cell Prep guide (document CG00053), which explains how to wash and resuspend cells.

10XG have an online database for publications, searchable by organism and tissue, useful to find protocols for specific tissues: <https://www.10xgenomics.com/resources/publications?query=&page=1>

CAUTION: obtaining nicely dissociated and alive cells is often the most important and tricky step of the whole experiment. Therefore, test and optimize the dissociation procedure beforehand, including any putative step where the cells will have to wait on ice before being processed by us (waiting at FACS, journey towards GECF...).

### Dissociation

- Perform at least one wash, as described in the 10XG cell prep guide. This wash can be performed with your culture medium. In case cell amount is too low, this can exceptionally be skipped (see cell prep guide).
- If thawing cells, follow either 10XG #CG000447 protocol or your preferred one, but include at least 1 wash step or “pseudo-wash” step (=thaw 1ml cells into 10ml medium → pellet → resuspend in chosen buffer).
- **Cell resuspension buffer:** cells should be brought either in PBS 0.04% BSA (1% for nuclei), or alternatively can be kept in their preferred medium (FCS – max 10% - is not a problem, as long as it was filtered to remove putative large particles). If you use BSA, use a molecular biology grade or similar. 10XG recommend:
  - Millipore/Sigma, BSA In DPBS, #A1595
  - ThermoFisher Scientific, UltraPure BSA, #AM2616
- When preparing cells/nuclei following protocols different from 10XG ones (e.g. Dead Cell Removal kit, Debris

Removal Solution, Miltenyi Anti-Nucleus MicroBeads) or for cells/nuclei coming from FACS, include a wash step to remove any additives and load the final resuspension in an approved buffer or media.

- Place cells on ice once prepared.
- Cells should be well dissociated from each other, with minimal amounts of doublets. If relevant carefully gate on FACS to avoid doublets.
- Absolutely avoid aggregates or clumps as they may clog the capillaries and lead to run failure. To avoid these, cells can be passed through a cell strainer, such as Flowmi cell strainer of 40um or 70um (we have both); or Miltenyi Biotec 30 um PreSeparation Filter #130041407. *Avoid using 30um filter if it comes with a risk of skewing your data.*
- Avoid cell debris as much as possible, as they are encapsulated along with cells and impact data negatively.
- Avoid using plastic petris or dishes for mincing/crushing tissues, as this was shown to create plastic debris that can clog the capillaries. Only use glass ones.
- If used, trypsin should be inactivated after use (serum, BSA...).
- EDTA is not OK here as it will impede reverse transcription. If you use EDTA to detach cells, quench it and do an additional wash. If really needed, limit it to 0.5mM, but still without guarantee of success. EGTA is better as it mostly chelates Ca<sup>++</sup> rather than Mg<sup>++</sup>.
- FACS buffer contains 1mM EDTA, so make sure the FACS buffer content in the final cells suspension is minimal.
- When obtaining clean single cells suspension is not possible, isolated nuclei can be run instead of cells:
  - 10XG offer a nuclei isolation kit (<https://www.10xgenomics.com/support/single-cell-atac/documentation/steps/sample-prep/chromium-nuclei-isolation-kit-sample-prep-user-guide>), we have it in house and can sell you aliquots for you to try.
  - 10XG also have several nuclei protocols on their website for specific cell types.
  - Nuclei tend to aggregate, especially if cytoplasm material is still attached or if some nuclear content extrudes. To reduce clumping: the recommended concentration of BSA in the wash and resuspension buffer is 1% instead of 0.04%. It can even be increased up to 2%. In parallel, PBS can be diluted to 0.7x (suggestion from 10XG tech support). Finally, nuclei concentration should be kept at a lower concentration range at all steps, and pelleting should be done at the lowest possible speed.
  - RNase inhibitor addition is mandatory at least in the wash and resuspension buffers, at 0.2 U/ul. See 10XG-recommended brands below. Also add some to FACS collection tube when relevant.
  - If the nuclei isolation protocol is not well established, a pilot is recommended, especially if it involves nuclei FACS sorting. Nuclei quality is indeed difficult to assess as one cannot check their “viability”, and they can be leaky. Alternatively, RNA can be extracted from nuclei and QCed (quantified/profiled).
- Notes
  - For dissociation of tricky tissues, refer to publications from the relevant field, contact 10XG support, look at users-developed protocols (<https://community.10xgenomics.com/t5/Customer-Developed-Protocols/ct-p/customer-protocols>) and/or refer to this 10XG webpage containing useful links: <https://kb.10xgenomics.com/hc/en-us/articles/218169563-How-do-I-dissociate-my-tissue-of-interest>, such as the Worthington Tissue Dissociation Guide. Using a GentleMACS instrument may help, as well as accutase, collagenase, pronase or psychrophilic proteases that work at 4° (Potter lab, PMID: 28851704).
  - Some people recommend DNase 1 treatment. For instance: <https://www.stemcell.com/how-to-reduce->

[cell-clumping-in-single-cell-suspensions-with-dnase.html](#) . Cells/tube walls must be especially well washed if used.

- Remember that expression profiles can change rapidly if dissociation is lengthy, in particular if not performed at cold temperatures. If you fear the dissociation procedure may alter the cell state, you can collect the dissociated cells, extract RNA and ask us to perform a bulk RNA-seq on these in parallel to the undissociated tissue of origin.
- Some tissues are specially RNases-rich, in particular pancreas, spleen and lungs (<https://kb.10xgenomics.com/hc/en-us/articles/4415486278669-Can-I-use-RNase-rich-tissue-samples-for-single-cell-gene-expression-or-Multiome-assays->). With such tissues, in particular if many cells are dying and release their content, it is recommended adding RNase inhibitor at least in the final resuspending buffer/medium (0.2-1U/ul depending on the protocol/application). 10XG recommend Sigma-Aldrich PN-3335399001. CAUTION: some inhibitors impact negatively performances:

**Other RNase inhibitors tested and shown to negatively impact data quality:**

- SUPERaseIn RNase Inhibitor
- Ribonucleoside Vanadyl Complex

*It is strongly recommended to use the RNase Inhibitor part number listed in the Multiome ATAC + Gene Expression nuclei isolation demonstrated protocols (Sigma Protector RNase inhibitor, PN-3335399001).*

Ambion™ RNase Inhibitor (AM2684 Thermo Fisher Scientific), RNaseOUT™ Recombinant Ribonuclease Inhibitor (10777019 Thermo Fisher Scientific) and, RiboLock RNase Inhibitor (Thermo Fisher, PN-EO0382) may be used as an alternative when Sigma Protector RNase inhibitor is not available. Note: there was limited testing of these inhibitors for Multiome ATAC + Gene Expression, so use with caution.

### Cell number and concentration

- Cells concentration should be 1'000-1'500 cells /ul. Concentration outside this range makes recovered cell number even less predictable than it is intrinsically. FACS often overestimates cells concentration, therefore take safety margin and ask the FACS facility to aim for a concentration a bit higher than that.
- We need > 50ul to perform the QC and the run itself in good condition for any number of targeted cells. If no more than 5'000 cells are targeted, 25ul could still be manageable if you cannot get more. If these values cannot be reached, contact us in advance.
- For concentrating or washing cells, 10XG recommend spinning at RT° at 300xg for 5min for small cells (5-10um), 250xg for 4min for medium cells (10-17um), and 150xg for 3 min for larger cells (17-25um). This ensures efficient pelleting while minimizing carryover of debris or aggregation of cells. When doing it the first time, keep and check also supernatant for absence of cells. If you have or suspect to have a sub-population composed of small cells, prefer small cells settings to avoid losing them.
- Before loading cells in the instrument, we'll double-check the cell concentration.
- Bring some extra resuspension buffer/medium to perform cells dilutions if needed.
- When cells tend to aggregate, keep their concentration on the lower side of the range.
- Automated counter such as the Countess are generally not accurate for very small cells.


### Viability

- The key word for the 10XG experiments is cell viability. It is very important to get good results downstream.

- Dead cells can generally be excluded bioinformatically as they contain a high percentage of mitochondrial RNAs. The percentage of dead cells that is considered acceptable depends on your experiment. If your cells come from a healthy suspension cell line, anything more than 10% dead cells is probably a bad sign, while if working from primary cells that underwent hours of dissection and sorting, 20% may be considered acceptable. We will move forward with any sample containing less than 20% dead cells, unless the user explicitly decides otherwise. When samples contain more dead cells, it is the user/PI decision to either move forward, perform a dead cells removal procedure (see below), or cancel the experiment (no fee).
- MACS-sorting in general leads to higher viability than FACS-sorting, so prefer the former if possible. If FACSing, prefer larger nozzle and lower pressure (discuss ahead with flow facility team). See 10XG CG000720 document for recommendations.
- If you sort the cells by FACS before loading, include a viability stain to get rid of dead cells if possible.
- Miltenyi offer a kit to remove dead cells if needed (“Dead Cell Removal kit”, 130-090-101), which works at least on mammalian cells (probably also on insect cells but to be confirmed).
- We'll check viability of the cells before loading them in the instrument.

## Number of reads/cell

- The choice of number of reads/cell impacts costs and therefore can influence how many cells you target, it is thus good to have a rough idea already starting. There are 3 factors that can influence this choice:
  - the size of the cells. For very small cells (5-10um), mRNA levels are low and sequencing saturates fast, thus 50-100k reads/cell capture most of the available information. In contrast big cells (>25um) contain a lot of mRNA can be sequenced up to 200k reads/cell and above and still capture new information.
  - the biological question: for clustering cells in groups of known cell types, and/or to delineate their broad transcriptional profile, 50k reads/cells is enough. If you want to zoom on specific genes, then the more reads the better (up to complete sequencing saturation, see above). If you also want to delineate new cell type with subtle differences, having more reads is also better.
- Number of genes detected directly depends on the cell size (=amount of mRNA molecules per cell). Sequencing depth chosen may thus be modulated according to cell size. See plot from 10XG (for v2 reagents):

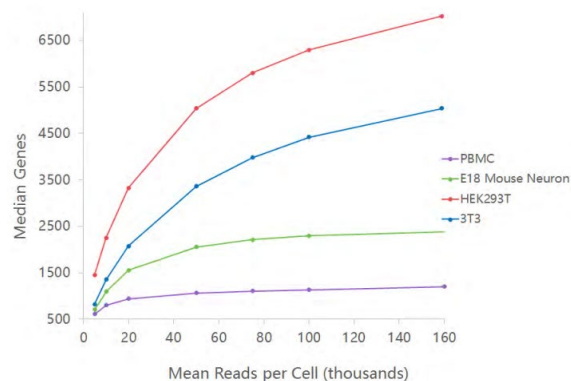

- 10XG recommend 50k reads/cells as a starting point. We rather go for 60-80k as a starting point.

## Diverse notes

- A general overview of 10XG experimental design: <https://www.10xgenomics.com/library/29cb3d>

- Special protocols allow to measure cell surface protein levels or monitor CRISPR perturbation screens along to the usual scRNA-seq. More information here: <https://www.10xgenomics.com/platforms/chromium>
- This method is intended for unfixed cells, but can be used with DSP/methanol-fixed cells as well (10XG protocol CG000776). For PFA-fixed cells, the FLEX method is possible (human and mouse only), see <https://www.10xgenomics.com/products/single-cell-gene-expression-flex> and inquire if interested.
- Samples must not contain toxic/carcinogenic chemicals. Contact us in advance if it is not possible.
- If TCR/BCR VDJ-seq is needed, tell us in advance as the whole method is different. See separated guidelines.
- If you have to prepare cells in different batches, make sure you avoid batch effects. Typically, avoid processing all controls in a batch and “treated” samples in a distinct batch.
- We can split the processing of cells in the instrument in several batches if this avoids some samples waiting on ice for very long times, and if this does not introduce batch effects.
- If you plan a big experiment, consider running just a single sample (pilot experiment) to assess quality of data and cell recovery rate.
- 10XG regularly update their reagents and workflows/pipelines. If you absolutely want us to use a specific version for comparing with a previously generated dataset, tell us ahead and we’ll discuss what can be done.
- Max cell diameter supported by 10XG is 30um, but the max theoretical limit is 65um (=size of the microcapillaries). These diameters are calculated for cells in suspension, not flattened on a petri dish. In case your cells are too big, use nuclei instead of cells (see above).
- When processing PBMC/blood cells, absence of left-over/contaminating erythrocytes is important to avoid losing reads to their ultra-abundant globin mRNA. In the same manner, if the isolation protocol includes an erythrocytes lysis step, it is critical to perform very careful washes afterwards to remove their floating mRNAs.
- Biosafety:
  - We process bsl2-level samples under a bsl2 hood. Please warn us in advance. Primary human cells must have been screened for absence of HIV, HBV and HCV infection (ideally also hCMV), or they will be considered bsl2. More generally, cells must not be infected with an infectious agent.
  - To determine the biosafety level of cell lines, this German website is used as a reference by the Swiss Biotechnology Office: <https://zag.bvl.bund.de/zelllinien/index.jsf?dswid=7026&dsrid=373>.
  - When nuclei are used, the isolation procedure may be sufficient to declassify your cells from bsl2 to bsl1, please investigate with us if relevant.
- Neutrophils are difficult to detect with this method (<https://kb.10xgenomics.com/s/article/360004024032-Can-I-process-neutrophils-or-other-granulocytes-using-10x-Single-Cell-applications>). Warn us in advance if this cell type is important for you, and rather go for the FLEX fixed-cells method (see our website).
- Neutrophils may not react well to long incubation on ice. They prefer room temperature during the sample prep. Other info about here: <https://kb.10xgenomics.com/hc/en-us/articles/360004024032-Can-I-process-neutrophils-or-other-granulocytes-using-10x-Single-Cell-applications->
- Samples can be multiplexed through an antibody-based system (TotalSeqB for 3’GE). See our separate guidelines for details.
- For the **bioinformatics** analysis, when an **ectopically expressed gene** must be added to the reference genome, give us its sequence in advance - including the 3’UTR since the reads will likely map there rather than in the coding sequence. If the transgene is expressed from a lentivector, the reads may lie in the 3’UTR that lies far away in the viral 3’LTR. Similarly, when tags/markers are inserted in 3’/C-ter of endogenous ORFs,

reads may lie within these added parts rather than within the endogenous gene. Finally, if the extra gene is already present in the genome, extra care is needed since multimapped reads are discarded by Cell Ranger.

- Droplets containing multiple beads occur at a mean frequency of ca 4%, and mRNAs from cells captured in these droplets are split into these multiple beads, hence will be detected as multiple cells, each with low UMIs. It is very likely that such “ghost cells” are filtered out during data processing, but to be on the safe side, new rare cell population composed of cells with low UMI counts should be validated by a different method.
- This is a polyA-based method, therefore it won't catch non-polyA RNAs.

## Versions log

- vA.01-vA.14: earlier releases. Indicated: that neutrophils are difficult to detect; the ghost cells issue; to pay special attention to erythrocytes when preparing PBMCs/blood cells; the minimum number of cells to bring; that when the transgene is expressed from a lentivector, the reads will be in the 3'LTR. Described specific procedures in case of bsl2 operations. Suggested to do bulk RNA-seq at the end of the dissociation procedure in case the procedure is harsh and may impact on cell states. Introduced CellPlex. Commented on negative impact of cell debris and clumps. Further considerations regarding biosafety. Clarified EDTA absence. Mentioned to use low pressure/large nozzle for FACS. Mentioned DNase 1 treatment for cells dissociation. Mentioned that only polyA RNAs are sequenced. Mentioned searchable publications database. Mentioned which BSA can be used to prepare the PBS 0.04% BSA. Added an important disclaimer section.
- vA.15: Clarified not to schedule an experiment without our prior agreement regarding a timeframe. Suggested the use of RNase inhibitors in some cases.
- vA.16: Mentioned: the submission form; to use glass recipients only when mincing/crushing tissues; Chromium Nuclei Isolation kit; Fixed cells (Flex method).
- vA.17: Recommendations to avoid nuclei clumps. And 1% BSA instead of 0.04 for nuclei. And RNase inhibitor recommended for nuclei.
- vA.18 (06.03.2024): Mentioned to use 1.5-2ml tubes to bring cells. Clarified that RNase inhibitor is mandatory for single nuclei, and the 0.2 U/ul concentration. Mentioned the new thawing protocol of 10XG, and the requirement for a wash when bringing thawed cells. Mentioned that neutrophils may prefer not to be kept on ice. Clarified that if an extra sequence must be added to the reference genome, if it contains an endogenous gene part, reads will be discarded due to multimapping.
- vA.19 (14.09.2024): Last guidelines version for v3.1 chemistry. Mentioned that an RNA QC is recommended when setting up complex nuclei isolation protocols. Mentioned the pilot in disclaimers. Mentioned that adding RNase inhibitor at extra steps may be required for lungs/spleen/pancreas.
- vA.20 (14.09.2024): First version specific for the new v4 chemistry, including several key changes (see first paragraph).
- vA.21 (19.11.2025): FACS buffer contains 1mM EDTA. Extra wash recommended when preparing cells/nuclei following protocols different from 10x Genomics ones (e.g. Dead Cell Removal kit, Debris Removal Solution or Miltenyi Anti-Nucleus MicroBeads). Added max 10% FBS. If the experiment is critical, and if users have extra cells, keep them on ice until the end of the run or fix them in DSP Methanol (CG000776).